

NUTRAKEM SURGICALS

Product Catalogue

www.nutrakemsurgical.com



About the Company

With 10+ years of industry experience, Nutrakem powers pharma with precision 4 from branded formulations to frontline hospital supply. We have a prolific presence in government procurement, forge strategic marketing alliances, and drive global export momentum. Our products carry the Nutrakem brand, crafted by trusted third-party manufacturers and marketed with impact. Our playbook? Compliance, positioning, and velocity 4 built for clarity, built for control.

Our vision is to become a global force in advancing healthcare by delivering world-class, affordable surgical and medical solutions that empower healing, inspire trust, and make quality care accessible to every corner of the world.

Our Memberships and Certifications


GOVERNMENT OF INDIA
MINISTRY OF CORPORATE AFFAIRS
Central Registration Centre
Certificate of Incorporation

(Pursuant to sub-section (2) of section 7 and sub-section (1) of section 8 of the Companies Act, 2013 (18 of 2013) and rule 18 of the Companies (Incorporation) Rules, 2014)

I hereby certify that NUTRAKEM INDIA PRIVATE LIMITED is incorporated on this FOURTH day of MARCH TWO THOUSAND TWENTY FOUR under the Companies Act, 2013 (18 of 2013) and that the company is Company limited by shares

The Corporate Identity Number of the company is **U46497MP2024PTC070135**
The Permanent Account Number (PAN) of the company is **AAJCN4661N***
The Tax Deduction and Collection Account Number (TAN) of the company is **BPLN06993A***

Given under my hand at Manesar this FOURTH day of MARCH TWO THOUSAND TWENTY FOUR

Signature Not Verified
Digitally signed by
DS MINISTRY OF CORPORATE
AFFAIRS, CIVIL MANESAR 1
Date: 2024.03.07 16:23:12 IST

Parikaj Srivastava
Assistant Registrar of Companies/ Deputy Registrar of Companies/ Registrar of Companies
For and on behalf of the Jurisdictional Registrar of Companies
Registrar of Companies
Central Registration Centre

Disclaimer: This certificate only evidences incorporation of the company on the basis of documents and declarations of the applicant(s). This certificate is neither a license nor permission to conduct business or solicit deposits or funds from public. Permission of sector regulator is necessary wherever required. Registration status and other details of the company can be verified on mca.gov.in
Mailing Address as per record available in Registrar of Companies office:
NUTRAKEM INDIA PRIVATE LIMITED
213, Sapphire Heights, Opp. C-21 Mall, A.B. Road, Vijay Nagar, Indore, Indore- 452010, Madhya Pradesh
*As issued by Income Tax Department


#startupindia

CERTIFICATE NO: DPP163429

GOVERNMENT OF INDIA
Ministry of Commerce & Industry
Department for Promotion of Industry and Internal Trade

CERTIFICATE OF RECOGNITION

This is to certify that NUTRAKEM INDIA PRIVATE LIMITED incorporated as a Private Limited Company on 04-03-2024, is recognized as a startup by the Department for Promotion of Industry and Internal Trade. The startup is working in "Healthcare & Lifesciences" Industry and "Pharmaceutical" sector as self-certified by them.

This certificate shall only be valid for the Entry up to Ten years from the date of its incorporation only if an turnover for any of the financial years has not exceeded ₹ 100 Cr.

21-04-2024 03-03-2024
DATE OF ISSUE VALID UPTO



Startup India


भारत सरकार
Government of India
सूक्ष्म, लघु एवं मध्यम उद्यम विभाग
Ministry of Micro, Small and Medium Enterprises

UDYAM REGISTRATION CERTIFICATE

UDYAM REGISTRATION NUMBER: UDYAM-MP-23-0140656

NAME OF ENTERPRISE: MUS NUTRAKEM INDIA PRIVATE LIMITED

TYPE OF ENTERPRISE *
S.No. Classification Year Enterprise Type Classification Date
1 2024-25 Micro 22/04/2024

MAJOR ACTIVITY: **TRADING**
(For availing benefits of Priority Sector Lending(PSL) ONLY)

SOCIAL CATEGORY OF ENTREPRENEUR: GENERAL

NAME OF UNIT(S)
S.No. Name of Unit(s)
1 NUTRAKEM INDIA PRIVATE LIMITED

OFFICIAL ADDRESS OF ENTERPRISE

Flat/Door/Block No.	SECOND FLOOR, 213	Name of Premises/ Building	SAPPHIRE HEIGHTS
Village/Town	INDORE	Block	SCHEME NO 54, OPP. C 21 MALL,
Road/Street/Lane	VIJAY NAGAR	City	INDORE
State	MADHYA PRADESH	District	INDORE , Pin 452010
Mobile	9103914720	Email:	nutrakemindia@gmail.com

DATE OF INCORPORATION / REGISTRATION OF ENTERPRISE: 04/03/2024

DATE OF COMMENCEMENT OF PRODUCTION/BUSINESS: 04/03/2024

NATIONAL INDUSTRY CLASSIFICATION CODE(S)

S.No.	NIC 2 Digit	NIC 4 Digit	NIC 5 Digit	Activity
1	46- Wholesale trade, except of motor vehicles and motorcycles	4649- Wholesale of other household goods	46497- Wholesale of pharmaceutical and medical goods	Trading

MSME Registration

Corporate Identification Number (CIN)


GOVERNMENT OF INDIA
Form GST REG-06
(See Rule 10(1))
Registration Certificate
Registration Number : 23AAJCN4661N122

1. Legal Name	NUTRAKEM INDIA PRIVATE LIMITED		
2. Trade Name, if any	NUTRAKEM INDIA PRIVATE LIMITED		
3. Additional trade names, if any			
4. Constitution of Business	Private Limited Company		
5. Address of Principal Place of Business	Floor No.: SECOND FLOOR Building No./Flat No.: 213 Name Of Premises/Building: SAPPHIRE HEIGHTS Road/Street: SCHEME NO 54, OPP. C 21 MALL, A.B. ROAD Locality/Sub Locality: Vijay Nagar City/Town/Village: Indore District: Indore State: Madhya Pradesh PIN Code: 452010		
6. Date of Liability			
7. Period of Validity	From	19/04/2024	To Not Applicable
8. Type of Registration	Regular		
9. Particulars of Approving	Madhya Pradesh		
Signature	Signature Not Verified Digitally signed by OFFICER IN CHARGE, GST REGISTRATION INDORE, MADHYA PRADESH Date: 2024.04.18 16:30:33 IST		
Name	KAILASH SULLIYA		
Designation	Assistant Commissioner of State Tax		
Jurisdictional Office	Indore - 11		
Date of issue of Certificate	19/04/2024		

Note: The registration certificate is required to be prominently displayed at all places of business in the State.

This is a system generated digitally signed Registration Certificate issued based on the approval of application granted on 19/04/2024 by the jurisdictional authority.

GST Registration


EXPORT PROMOTION COUNCIL FOR MEDICAL DEVICES
Export Promotion Council for Medical Devices
205, Second Floor, Commercial Complex, Omega-1 Greater Noida, Noida,
GAUTAM BUDDHA NAGAR, UTTAR PRADESH, 201110
Email: membership.officer@epcmid.in, Telephone: +91 9011817661 Website: <http://epcmid.in>

REGISTRATION - CUM - MEMBERSHIP CERTIFICATE
(Valid under the provisions of Foreign Trade Policy: Government of India)

Name of the exporter	NUTRAKEM INDIA PRIVATE LIMITED		
IEC Number	AAJCN4661N	PAN	AAJCN4661N
Date of Birth / Incorporation	04/03/2024		
Address of the Registered Head Office	213, Sapphire Heights, Opp. C-21 Mall, A.B. Road Vijay Nagar, Indore - 452010 (Madhya Pradesh) Indore, 452010, INDORE, MADHYA PRADESH		
Nature of Concern	Private Limited		
Description of goods/ services for which registered	All medical devices and products		
Registration Number	ECMC/EPCMD/0013/2023-2026		
Issued against File Number	ECMC/APPLY/EPCMD/0001/068AM26		
Registered as	Merchant Exporter		
Name of the Proprietor/ Partner/Director/Shareholder	Akshita Pandey-Anil Kumar Pandey		

This certificate is issued as per the details of our records and is subject to the conditions laid down in the relevant scheme of registration of this council.

Issue Date: 18/09/2023 Valid upto: 31/03/2026

Note: Refer DGFT Website (Services - View IEC Related Details - View App IEC) for IEC Branch and Director details.

Validity unknown

Digitally signed by
Date: 23-Sep-2023 16:30:33 IST
Reason: operation of system-generated in Location: Export Promotion Council for Medical Devices (EPCMD)

This document is digitally signed by membership.officer@epcmid.in (Praveen Kumar Mittal), on 22/09/2023.

Export Promotion Council for Medical Devices- Membership Certificate


GOVERNMENT OF INDIA / MSME विभाग
Ministry of Commerce and Industry / सूक्ष्म एवं मध्यम उद्यम विभाग
Directorate General of Foreign Trade / विदेश व्यापार निदेशकालय
Office of the Joint Director General of Foreign Trade, Indore अंशक निदेशक विदेश व्यापार का कार्यालय, इंदौर
BSNL Building, One CTO Compound, Residency Area, INDORE, MADHYA PRADESH, 452001 /
भारतव्यापार निदेशक, एक सीटीओ कंपाउंड, रेसिडेंसी एरिया, इंदौर, इंदौर, मध्य प्रदेश, 452001

This is to certify that NUTRAKEM INDIA PRIVATE LIMITED is issued an Importer/Exporter Code (IEC) AAJCN4661N with details as follows:-

IEC	AAJCN4661N
संयुक्त निर्यात कोड (IEC) / PAN	AAJCN4661N
फर्म का नाम / Firm Name	NUTRAKEM INDIA PRIVATE LIMITED
व्यवसाय की प्रकृति / Nature of Concern	Private Limited
जारी करने की तिथि / Date of Issue	18/09/2023
पता / Registered Address	213, Sapphire Heights, Opp. C-21 Mall, A.B. Road, Vijay Nagar, Indore - 452010 (Madhya Pradesh), Indore, INDORE, MADHYA PRADESH, 452010
सहस्रनामा / Name of the Signatory	Akshita Pandey
निदेशक / Partner Details	Refer online at https://dgtf.gov.in or scan the QR Code
संयुक्त निर्यात कोड / Branch Details	Refer online at https://dgtf.gov.in or scan the QR Code

Last Modified : 18/09/2023
File Number : IENR/EC/PPLY/0000/992AM26

Note : This is a system-generated certificate. Authenticity / Updated details of the IEC can be checked at official DGFT website <https://dgtf.gov.in> by entering the IEC and Firm Name under Services > View App IEC Details. You can also authenticate the certificate by scanning the QR code.

Import Export Code

Why Choose Us?

Quality Standards

- Products made in compliance international standards
- WHO-GMP, ISO 13485 & 9001 Compliance
- CE Certified, MDR Approved

Manufacturing Excellence

- Third party manufacturing facility available
- Works on Standard SOP
- Made with International Standards

Testing & Delivery

- Inbuilt Sterility Testing Lab
- Inbuilt Toxicity Testing Lab
- Inbuilt Leakage/Sealing Testing Lab
- On time and fast Delivery
- Competitive rates

Our products have the following certifications



Surgical Sutures

- Hightensile strength
- Synthetic
- Hermetically sealed packing
- Needle type: Taper Point, Straight, Cutting Edge, Blunt Point, Round Bodied



Absorbable Sutures



Polyglycolic Acid Suture

- Best use for general soft tissue closing and/or ligation
- Used in general surgery, gastrointestinal surgery, ophthalmic surgery, plastic surgery, skin closure, obstetrics, gynecology urology, and orthopedics



Polyglactin 910 Suture

- Synthetic absorbable sterile surgical suture composed of a copolymer
- Material: 90% Glycolide and 10% L-lactide.
- Material break down by hydrolytic action results in total absorption in approx 56 - 70 Days
- Approx 50% Tensile Strength remains After 7 Days
- Tensile strength: 75% for 14 Days and 50% till 21 Days

Poliglecaprone 25 suture



- Monofilament synthetic absorbable surgical suture
- Excellent tensile strength and smooth in handling
- Material: Poliglecaprone 25, A Copolymer Made Of 75% Glycolide And 25% Caprolactone
- Poliglecaprone suture is ideal for fast-healing tissues, ensuring reliable and effective wound closure
- Nonantigenic, non pyrogenic and elicits only a slight tissue reaction during absorption
- Hydrolytic action results in total absorption in approximately 90 to 120 days.
- Tensile Strength: Approx 60% Of tensile strength remains After 7 Days and 50% after 14 days.

Monofilament Polydioxanone Suture (PDS)



- Material: Prepared from the polyester, poly(p-dioxanone)
- Minimizes variability of tensile strength retention and absorption rate providing wound support through an extended healing period
- Hydrolytic Action results In total absorption in approx 180 - 210 Days
- Tensile Strength: Approx 50% tensile strength remains After 42 Days ; 50% Tensile Strength Is Retained For 14 Days

Chronic Catgut Suture

- Absorbable, sterile surgical suture composed of purified connective tissue (mostly collagen)
- Derived from the sub mucosal fibrous layer of healthy sheep intestines
- Provides greater resistance to absorption
- Suitable for use in general, surgical procedures, in soft tissues repairing procedures, or in gastrointestinal surgery
- Tensile Strength Is Retained For 14 - 21 Days



Plain Catgut Suture

- Fast absorbing surgical gut suture
- Material: Made from strand of collagenous material prepared from the sub mucosal layers of the small intestine of healthy sheep
- is sterile and elicit slight tissue reaction during absorption
- Absorption by phagocytosis: Complete mass absorption In approx 65 Days.
- Tensile Strength: Tensile Strength Is Retained For 7 - 14 Days. Approximately 50% Of Tensile Strength Remains After 7 Days.
- Suitable for surgical procedures for tissue anastomosis and ligature and have low tissue reactivity.



- Absorbs more quickly when used in infected tissues or tissues with increased proteolytic enzyme levels (Stomach, Cervix, Vagina)


Non-Absorbable Sutures



- Preferable for critical surgeries - Offers high tensile strength, minimal tissue reaction, and excellent knot security.

Monofilament Polypropylene Suture

- Sterile surgical non-absorbable suture
- Material: Composed of an isotactic crystalline stereoisomer of polypropylene, a synthetic linear polyolefin.
- Pigmented blue to enhance visibility
- Used in general soft tissue approximation and/or ligation as well as cardiovascular, ophthalmic, and neurological procedures.
- Excellent tensile strength and knot tying with minimum tissue reaction

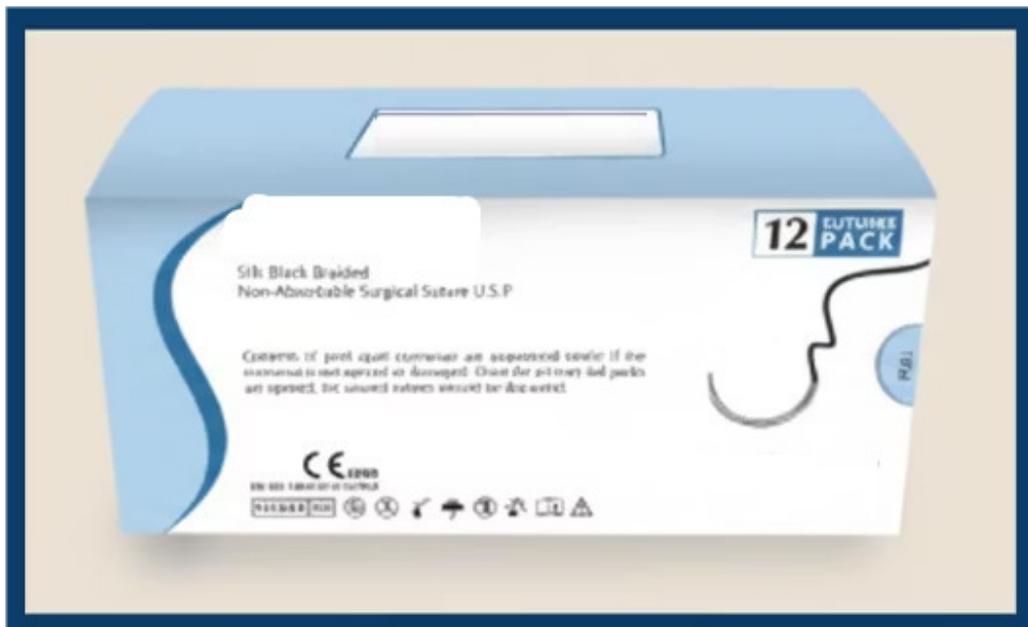

- Absorption: An early inflammatory reaction in tissue is followed gradual encapsulation of suture by fibrous connective tissues

Black Braided Silk Suture

- Sterile surgical non absorbable suture
- Material: Composed of an organic protein called fibroin derived from domesticated species Bombyx mori (B. More) of the family Bombycidae.
- It is obtained from the braided fibres from the cocoon of the silkworm coated with wax silicone
- Tensile Strength: Excellent tensile strength and knot tying with minimum tissue reaction
- Uses in general soft tissue approximation and/or ligations as well as cardiovascular, ophthalmic and neurological procedures.

Virgin Silk Suture

- Sterile surgical non absorbable suture
- Material: Made from filaments of silk fibre held together by natural gum secreted by silk worm i.e., twisted fibres from the cocoon of the silkworm
- Absorption: Silk Suture Elicits an Initial Inflammatory Reaction In Tissues, Which Is Followed By Gradual Encapsulation Of The Suture By Fibrous Connective Tissues.
- Tensile Strength: Excellent tensile strength and knot tying, minimum tissue reaction


- Used in in general soft tissue approximation and/or ligation, including as well as cardiovascular, ophthalmic and neurological procedures.



Braided Polyester Suture

- Sterile surgical non absorbable suture
- Material: Composed of Polyester - Polyethylene terephthalate coated with silicone
- Highly adherent coating, relatively nonreactive and acts as a lubricant to mechanically improve physical properties of uncoated suture by improving handling qualities
- Braided for optimal handling
- Green dyed for enhanced visibility
- Tensile Strength: Excellent tensile strength and knot tying, minimum tissue reaction


- Absorption: Non-absorbable, gradually encapsulated by connective tissue

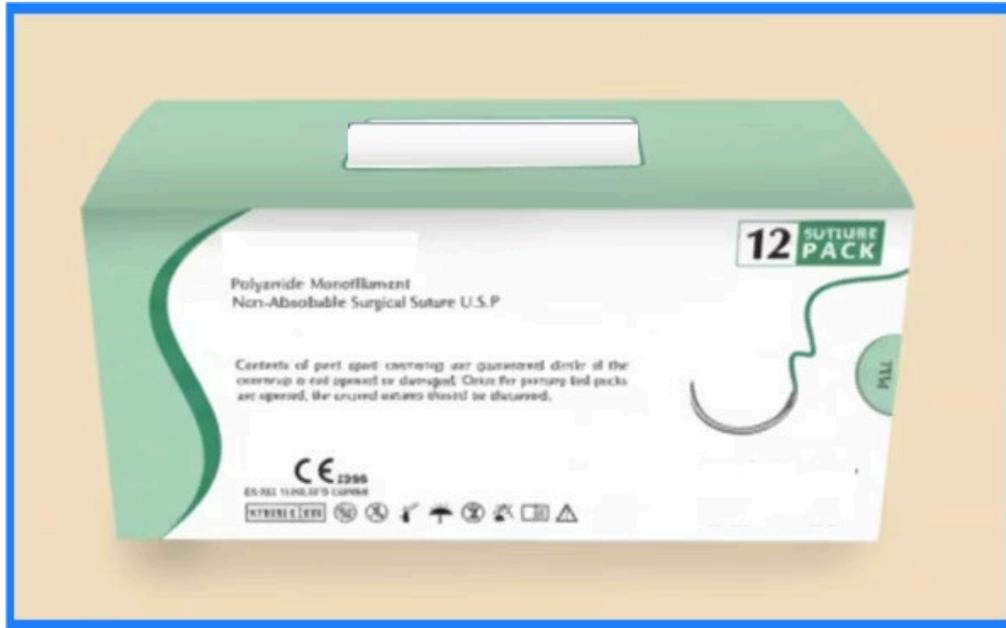

- Commonly referred to as nylon suture its use is widespread to its high tensile strength, non-absorbability and elasticity

Monofilament Polyamide Suture

- Material: Composed of the long-chain aliphatic polymers Nylon 6
- Dyed black / blue / green to enhance visibility in tissue
- Absorption: Non-absorbable, gradually encapsulated by connective tissue.
- Thread mass diminishes, approx 10% a year by rupture of chemical links through hydrolytic action
- Tensile Strength: Excellent tensile strength and knot tying, minimum tissue reaction
- Mostly used in ophthalmic and neurological procedures
- Sterile surgical non absorbable suture

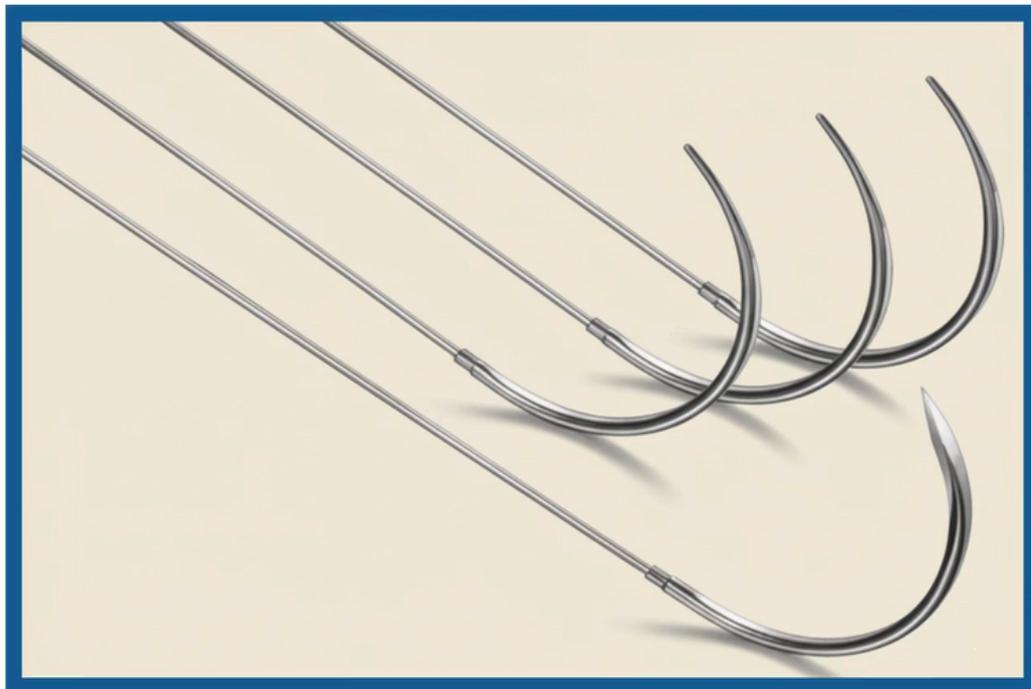

- Monofilament and multifilament suture

Stainless Steel Suture

- Sterile surgical non absorbable suture
- Material: Composed of 316L stainless steel conforming to ASTM Standard; F138 Grade 2 Stainless Steel bar; and wire for surgical implants
- Meets United States pharmacopeia (USP) Standards
- Used surgeries requiring durable wound support such as orthopedic, cardiovascular, and general surgeries
- Tensile Strength: Excellent tensile strength and knot tying, minimum tissue reaction
- Absorption: Non-absorbable

Mesh & Bone wax



Surgical Bone Wax

- Mix of beeswax, paraffin, and isopropyl palmitate, a wax-softening agent.
- Has a waxy odour and is opaque
- Used for control bleeding from bone surfaces during surgeries
- Available in sterile individual foil envelopes, of 2.5 g and packed in individually sealed overwrap, sterilized by Gamma radiation.
- Absorption: Minimally Resorbable
- Uses: Orthopaedic Surgery, Neurosurgery, Thoracic Surgery, Traumatology, Maxillofacial Surgery, Thoracic Surgery, Dental, Oral and Maxillary-Face Surgery
- Its non-absorbable and non reactive biochemically


- Ensures hemostasis of bone by becoming a mechanical (tamponade) barrier.



Monofilament Polypropylene Mesh

- Light Weight Construction ensures less foreign mass in the body & flexible scar tissue
- Non Sterile

- Material: Made from knitted filament of extruded polypropylene
- This bi-directional elastic property allows adaptation to various stresses encountered in body
- Available in various thicknesses, dimensions, knitting patterns, shapes, sling, pore sizes and pre-cut sizes
- Type and shape of mesh available as per requirement
- Product features: Smooth, strong, thin and flexible, macro porous for high visibility and colonization
- Used for surgical procedures, pelvic floor reconstruction, tissue carrier, vaginal prolapse suspension, hernia procedures, urinary incontinence sling
- Rounded at both ends, smooth, flat and thin
- Packed in loose

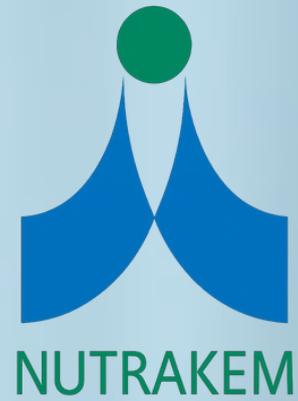

NUTRAKEM SURGICALS

LET'S WORK TOGETHER!

213 SapphireHeights, OppC-21 Mall,

Vijay Nagar, Indore, M.P. (India) - 452010

Mobile: +91-8103914720, +91-9687052275 | Email: info@nutrakemsurgicals.com

www.nutrakemsurgicals.com